

Administrator Manual

1. To access the system, you must:

- Enter your username.
- Enter your password.
- Click *Login*.

If your credentials are correct, you will be redirected to the main screen where you can upload radiographs.



Enter to your account

User

Password


Log In

2. On the main screen, you will find the **Records** of the classifications.

- Click the button **Export** to download general information about the records.
- Click the button **Download PDF** to download a specific record with its heatmap
- Click the button **Delete** to delete all the records.


Classifications Graph

- A bar chart displays the number of classifications daily or weekly.
- This helps analyze system activity and model usage over time.

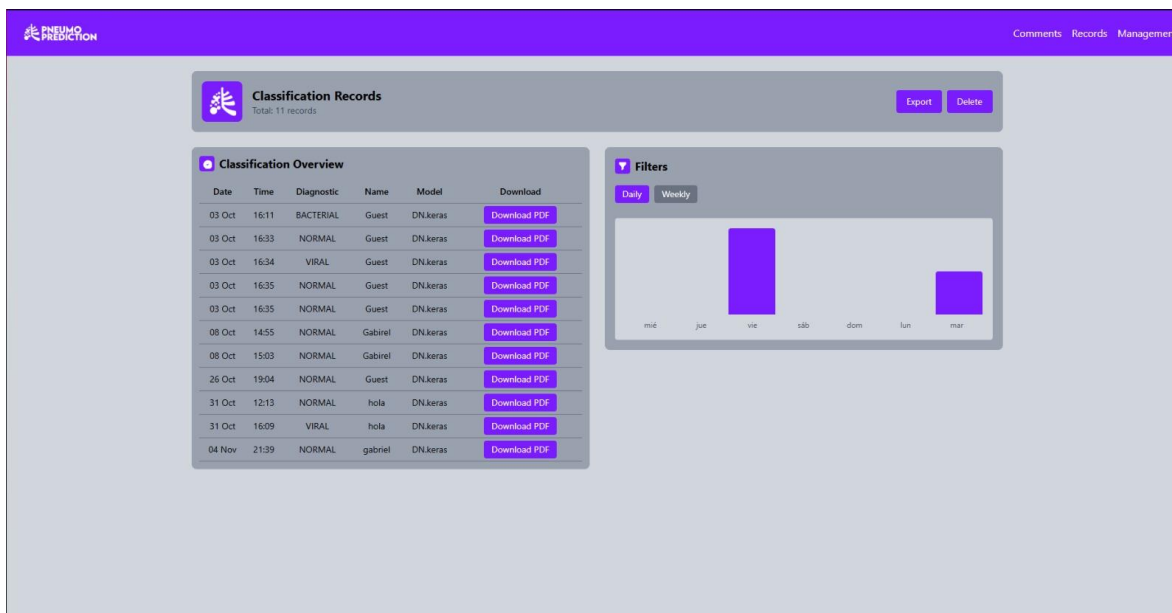

3. In the top navigation bar, you will find the Comments section, here you can:

- Read the comments that were created in the support section by the users.
- Input a specific word in the search bar to search for any comment with it.
- Delete the comment if it is needed.

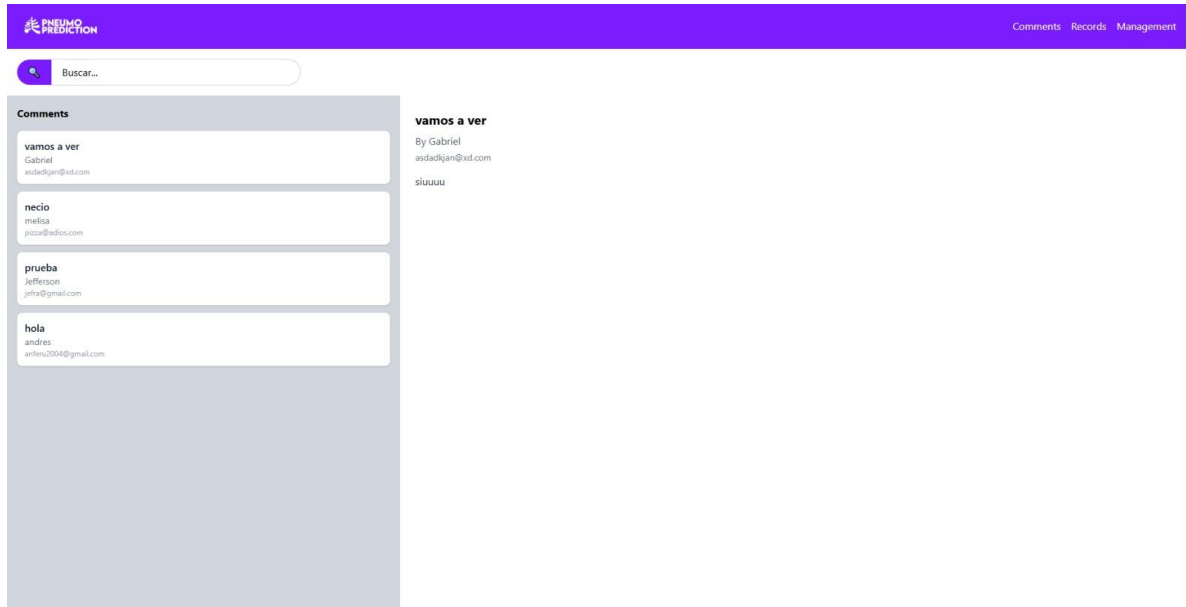

4. In the top navigation bar, you will find the Management section, here you can
- Create, delete or update any medic account.
 - Change the model that is being used in the classification.
 - Log out of your account safely using the Logout option.

The screenshot displays the 'PNEUMO PREDICTION' interface with a top navigation bar containing 'Comments', 'Records', and 'Management'. The 'Management' section is active, showing two panels: 'Account Management' and 'Model Configuration'.

Account Management
Manage your user accounts.

Create / Edit account

Name

Username

Password

Save

Modify account

Select user

Edit **Delete**

Model Configuration
Adjust the classification model.

Select model

Upload

Log out
Exit your account

Logout